



**TD-100™**  
**Transesophageal (TEE) Ultrasound Probe**  
**Disinfector**

**SERVICE MANUAL**

Manufactured by:



3300 Lake Woodard Dr.  
Raleigh, NC 27604 USA

[www.csmedicalllc.com](http://www.csmedicalllc.com)

Sales: [sales@csmedicalllc.com](mailto:sales@csmedicalllc.com)

Customer service: [service@csmedicalllc.com](mailto:service@csmedicalllc.com)

Phone: (877) 255-9472

Phone: (919) 255-9472

Fax: (919) 212-8050

TD-100™ is a registered trademark of CS Medical LLC.

# TD-100™ Probe Disinfector Operators Manual

## Table of Contents

<b>Section</b>		<b>Pages</b>
<b>1 TD-100™ System</b>		
1.1	Signal Words/Definitions	4
1.2	Illustration of the TD-100™ Probe Disinfector	6
1.3	System Overview of TD-100™ Probe Disinfector	7
1.4	Site Specifications	9
	-Electrical Requirements	
	-Water Requirements	
	-Environmental Requirements	
	-Classification	
<b>2</b>	<b>TD-5™ High Level Disinfectant</b>	
2.1	Vapor Control	11
2.2	Emergency Shut Down Procedures	12
<b>3</b>	<b>Operating Instructions</b>	
3.1	Aborted Cycle	13
<b>4</b>	<b>Quality Assurance</b>	
4.1	Valid Cycle Printout	14
4.2	Aborted/Invalid Cycle Printout	15
4.3	Troubleshooting	16
4.4	Water Supply Problem	
4.5	Diagnostic Failure	17
4.6	Trouble Shooting	18

<b>Section</b>		<b>Pages</b>
<b>5</b>	<b>Maintenance</b>	
	5.1 Daily Checks and Daily Cleaning	21
	5.2 Changing the Printer Paper	22
	5.3 Changing the Air Filter	23
	5.4 Setting or Changing the Time and Date	24
<b>6</b>	<b>Warning and Caution Statements</b>	
	6.1 Warning Statements	25
<b>7</b>	<b>Spare Parts</b>	
	7.1 Spare Parts List	26
<b>8</b>	<b>Contact Information</b>	
	8.1 Contact Information	27

# SECTION 1: TD-100™ System

## Section 1.1: Signal Words/Definitions



**WARNING!** - INDICATES A POTENTIALLY HAZARDOUS SITUATION, WHICH IF NOT AVOIDED, COULD RESULT IN DEATH OR SERIOUS INJURY.



**CAUTION!** - INDICATES A POTENTIALLY HAZARDOUS SITUATION, WHICH IF NOT AVOIDED, MAY RESULT IN MINOR OR MODERATE INJURY. IT MAY ALSO BE USED TO ALERT AGAINST UNSAFE PRACTICES OR POTENTIAL EQUIPMENT DAMAGE.

## Definitions

**Purge:** Purge is when the unit is drained prior to starting a run, during normal operation or when a run is interrupted.

**TD-5™:** TD-5™ is the bottle which contains the disinfectant.

## Symbols



Symbol for Alternating Current



## **WARNING!**



**THE TD-100™ IS DESIGNED FOR USE WITH TEE ULTRASOUND PROBES ONLY.**



**THESE INSTRUCTIONS MUST BE PLACED IN A VISIBLE LOCATION.**



**DO NOT USE THIS SYSTEM UNTIL IT HAS BEEN PROPERLY INSTALLED.**



**THE TD-100™ SYSTEM MAY NOT BE USED IN EXPLOSION ENDANGERED AREAS.**



**CONNECT THE SYSTEM TO A HOSPITAL GRADE, GFI RECEPTACLE ONLY. GFI OUTLET MUST BE TESTED ON A REGULAR BASIS PER THE GFI MANUFACTURERS' RECOMMENDATION.**



**DO NOT ATTEMPT TO ACCESS NON-USER SERVICEABLE COMPONENTS OF THE TD-100™ SYSTEM.**



**DO NOT OPERATE THIS UNIT IF THE POWER CORD APPEARS DAMAGED OR FRAYED IN ANY WAY.**



## **CAUTION!**



**ONLY ULTRASOUND PROBES LISTED IN THIS MANUAL MAY BE USED WITH THIS SYSTEM - SEE LIST IN SECTION 9.2.**



**DO NOT OPERATE THIS UNIT IF DISPLAY IS BLANK OR HARD TO READ. IF DISPLAY IS HARD TO READ, CALL CUSTOMER SERVICE.**

# Illustration of the TD-100™ Probe Disinfecter

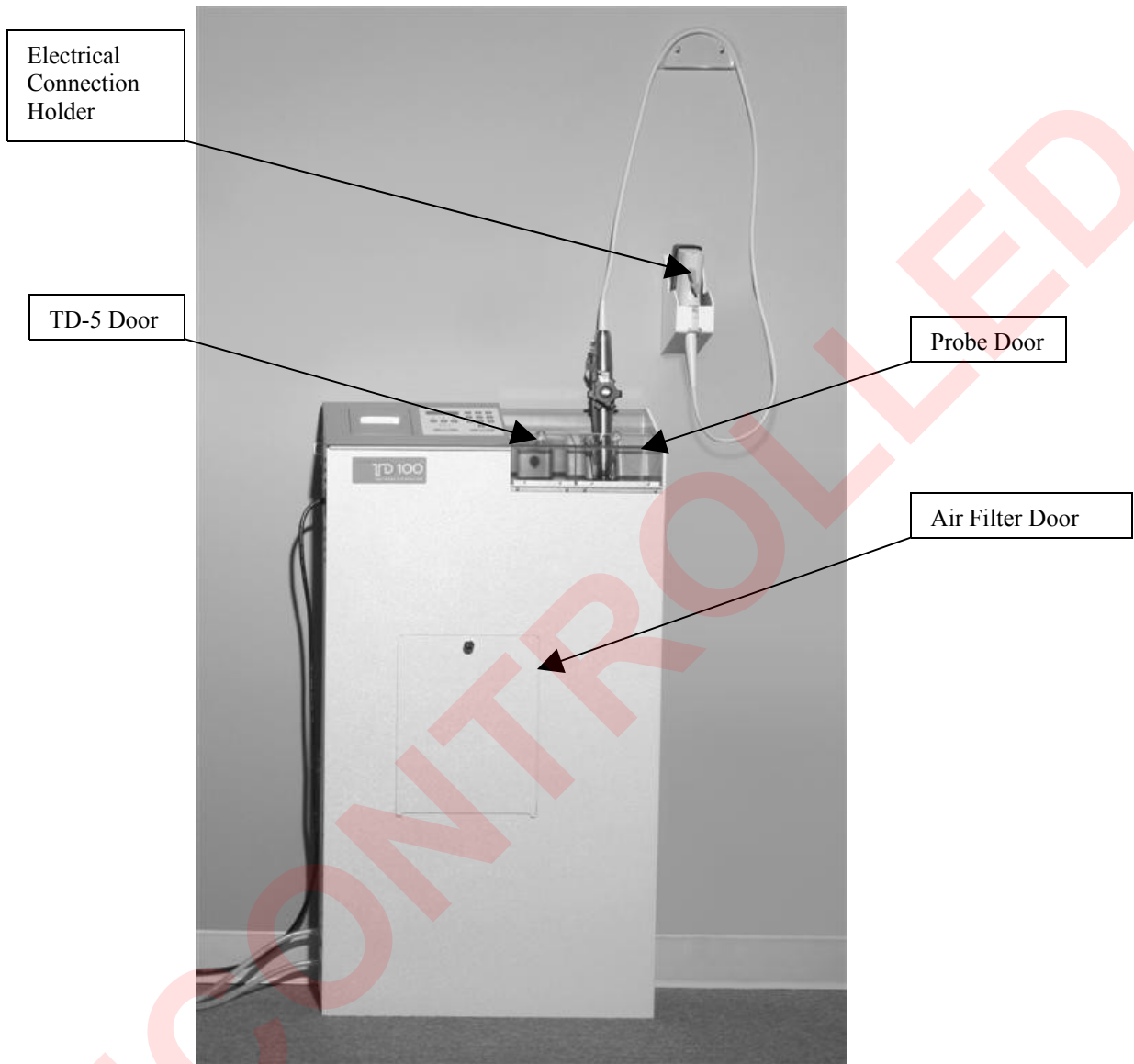


Fig. 1 Typical installation of the TD-100™.

### 1.3 System Overview of TD-100™ Probe Disinfectant

The TD-100™ is an automated disinfectant that is designed to provide high-level disinfection of transesophageal (TEE) ultrasound probes. A TEE ultrasound probe is considered a semi-critical device that enters the body through the esophagus. The TD-100™ rinses the disinfected TEE probe using high quality potable water.

**Note:**

When using potable water for rinsing, the user should be aware of the increased risk of recontaminating the TEE probe with *Pseudomonas* and atypical (fast growing) *Mycobacteria* often present in potable water supplies. Additionally, *Mycobacteria* are highly resistant to drying; therefore, rapid drying will avoid possible colonization but may not result in a device free from atypical *Mycobacteria*.

Water treatment systems, such as water softeners or deionizers, may add microorganisms to the treated water to the extent that microbial content of the water at the point of use could exceed that of the pretreated drinking water. To ensure proper water quality, adherence to maintenance of the water treatment system is recommended.

The use of a bacterial retentive (0.2 micron) filter system may eliminate or greatly reduce the amount of waterborne bacteria from the potable water source. Contact the manufacturer of the filter system for instructions on preventive maintenance to avoid colonization or formation of biofilms in the filter.

Endoscopes with their intricate channels can provide ideal situations for rapid colonization of bacteria. Although TEE probes do not have any internal channels, there is still the potential for colonization of waterborne bacteria. Waterborne bacteria are highly resistant to drying; rapid drying will avoid possible colonization. A final drying step is recommended following the TEE probe manufacturer's recommendations. In the absence of a recommended drying procedure the following is recommended: A final wipe down using 70% isopropyl alcohol solution is useful to speed the drying process and reduce the numbers of any organism present as a result of rinsing with potable water.

Although these bacteria are not normally pathogenic in patients with health immune systems, AIDS patients or other immuno-compromised individuals may be placed at high risk of infection by these opportunistic microorganisms.

The system uses TD-5™ disinfectant, which is designed only for use with the TD-100™. The disinfectant is packaged in 500ml bottles. One bottle of disinfectant is used for each process cycle. The disinfectant bottles are not designed to be re-used in the system. It is a violation of Federal law to use this product in a manner inconsistent with its labeling.

During each disinfection cycle, the machine itself is disinfected with the TD-5™ disinfectant. All wetted surfaces and components within the TD-100™ are exposed to the same disinfection conditions as TEE probes. After being disinfected these same surfaces, components and TEE probes are rinsed free of any residual TD-5™ disinfect, completing the disinfection cycle.

The TD-100™ is microprocessor controlled and will automatically disinfect the probe once the cycle begins.



**WARNING! TEE PROBES MUST BE PRE-CLEANED PRIOR TO USE IN THE SYSTEM. SEE PRE-CLEAN INSTRUCTIONS IN SECTION 3 OF THIS MANUAL.**

The TD-100™ will take approximately 17 minutes to complete the disinfection and rinse cycles, after which the TEE probe is removed from the machine.

At the end of a cycle, a printout will allow the user to verify that the machine has completed a valid cycle. If an error occurs during the process the run is invalidated and an error report prints.

At the end of a valid cycle, the TEE probe must be removed from the system and dried according to the probe manufacturer's recommendations.

The unit is placed on the floor against a wall. The installation kit includes four wall brackets for support. It requires an electrical connection to a 120V outlet, water and drain connections.



UNCONTROLLED

## 1.4 Site Requirements and Specifications

### Actual Size

Width	24"	[61cm]
Depth	8"	[20.32cm]
Height	44"	[112cm]
Weight	70 lbs	[32kg]

### Electrical Requirements

Service	120 V~ (VAC), 20 Amp, 60 HZ- dedicated circuit terminated in 20 Amp hospital grade GFCI double wall receptacle
Ratings	120 V~, 60 Hz, 8A
Fusing	Fast Acting, 10A, 125V, 5mm X 20 mm Cooper Bussman GMA-10A (or equivalent)

### Water Supply Requirements

Delivery pressure	Regulated to 20-25 PSI Using a Watts 263A with a pressure gauge and female 1/4" NPT outlet (or equivalent)
Flow rate	1gpm at 20 PSI Minimum
Temperature	Cold water supply: 50 - 80 degrees F
Shut-off	A shut-off valve is recommended
Connection	10 feet of 1/4" ID Tubing (provided in installation kit)

Quality	High quality potable water; a high quality potable water is one that meets Federal Clean Water Standards at the point of use.
Water usage	Approx one gallon per processing cycle

### **Drain Requirements**

Drain	1½” Drain pipe no more than 18” above the floor. A special drain connection fitting is provided in the installation kit.
-------	--

### **Environmental Requirements**

- Indoor use
- Altitude up to 2000m (6562 ft)
- Maximum relative humidity 80% for temperatures up to 31°C, decreasing linearly to 50% at 40°C.
- Mains supply fluctuations up to ±10% of the nominal voltage
- Category II installation
- Pollution degree 2 environment
- Room Temperature 60 – 90°F
- Room Humidity 20 - 80% RH non-condensing

### **Classification**

Protective against electric shock: Protective Class 1, ordinary protection.

## SECTION 2: TD-5™ High Level Disinfectant

### 2.1 Vapor Control

The American Congress of Industrial Hygienists (ACGIH) recommends a glutaraldehyde ceiling exposure limit of 0.05 ppm.

The TD-100™ is fitted with an air filtration system that minimizes exposure to glutaraldehyde fumes and vapors. During factory testing, normal exposure levels were found to be below the 0.05 ppm exposure limit. Correct use of this equipment per the guidelines and instructions in this Operators Manual will ensure that your exposure to glutaraldehyde vapors is minimized.



**CAUTION! SHOULD YOU SMELL GLUTARALDEHYDE VAPORS DURING NORMAL OPERATION, IMMEDIATELY CHECK EXHAUST FAN OPERATION BY PLACING HAND OVER OUTLET ON LEFT SIDE OF MACHINE AND FEEL FOR AIR FLOW. THEN SHUT DOWN THE SYSTEM AND FOLLOW THE EMERGENCY PURGE PROCEDURE.**

#### **EMERGENCY PURGE PROCEDURE**

- 1) Turn off main power switch for 10 seconds and turn back on.
- 2) System will purge any liquid in reservoir.
- 3) Once reservoir is purged, open the probe door and slowly pour in 16 oz of clean water.
- 4) Allow system to purge water and repeat step 3 with another 16oz of clean water.
- 5) Allow system to purge and shut off main power switch.
- 6) Remove TD-5 bottle, replace cap and dispose of.
- 7) Remove probe and dry according to probe manufactures' recommendation.



**WARNING! PROBE IS NOT DISINFECTED AND SHOULD NOT BE USED UNTIL PROPERLY DISINFECTED.**

## **2.2 Emergency Shut Down Procedure**

- 1) Turn the main power switch to the off position.
- 2) Turn off water.

### **AIR FILTER OR FAN FAILURE**

If you smell glutaraldehyde vapors and are following correct operating procedures per this manual; then either the fan is not working or the air filter has reached saturation.

- 1) Fan failure: Check airflow with power on by placing hand over outlet on left side of machine. If no air, call customer service.
- 2) Filter saturation: Change Filter (refer to Section 6.3)

## SECTION 3: Operating Instructions

### 3.1 Aborted Cycle

If during the run the unit detects an error, an abort message will display. The unit will display:

```
RUN ABORTED
YES REPRINT NO
```

The unit will print an invalid run report and request if you want a reprint of this report. Should you run out of paper, replace paper and press “A” to reprint the aborted run. The unit will always print “**WARNING PROBE NOT DISINFECTED**” whenever an invalid run occurs for any reason.

```
TD-100
RUN ABORTED
*****
SYSTEM FAIL
WARNING
PROBE NOT DISINFECTED
*****

01/11/06
11:15:51A

PROBE #: 0
OPERATOR #: 0
TD-5 LDT #: 000

SYSTEM TEST: FAILED
REMARKS:
_____
_____
_____
_____
_____
_____
_____
_____

DIAG V1018:
01/11/06
11:15:51A
MODE: 021 000 004 000
A IN: 253 011 171 193 223
INIT: 253 011 000 171 194
TIME: 001 101 000 000

HIST: 000 045 041
M0: 001 020
M3: 004 000 000
M4: 006 004 008
M5: 003 000 000
M6: 024 000 000 002 006
M7: 000 000
```



**WARNING! THE TEE PROBE IS NOT DISINFECTED WHENEVER AN INVALID RUN OCCURS. THE PROBE IS CONSIDERED CONTAMINATED UNTIL A VALID RUN WITH A VALID RUN PRINTOUT OCCURS.**

The printout on the left indicates an invalid run has occurred. The system prints the date and time of the invalid run and always lists the Diagnostics for this run and the history of previous runs. The diagnostic information is used to assist in troubleshooting the problem. When you call for assistance, please have a copy of your invalid run printout.



**CAUTION! ALWAYS SAVE THIS PRINTOUT FOR LATER REFERENCE.**




**CAUTION! ALWAYS RESET THE TD-100™ AFTER AN ABORTED CYCLE.**

To reset the machine, turn the power off and then turn the power on. This ensures that the TD-100™ purges any liquid left in the system after an aborted run. If after resetting the machine it is not possible to obtain a valid cycle, purge the machine using the EMERGENCY PURGE PROCEDURE (Section 2.4) and call customer service.

## SECTION 4: Quality Assurance

### 4.1 Valid Cycle Printout

This picture shows a valid cycle printout.

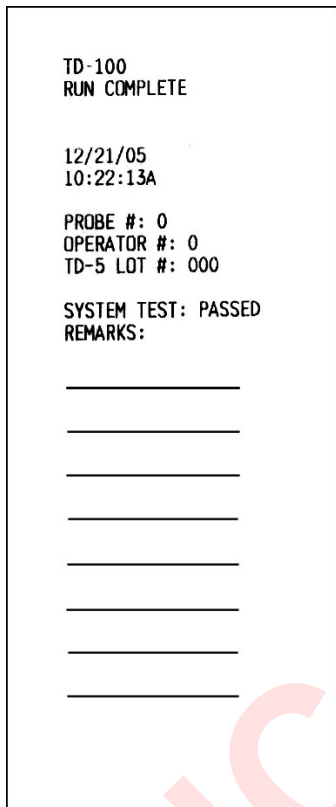

```
TD-100
RUN COMPLETE

12/21/05
10:22:13A

PROBE #: 0
OPERATOR #: 0
TD-5 LOT #: 000

SYSTEM TEST: PASSED
REMARKS:
_____
_____
_____
_____
_____
_____
_____
_____
_____
_____
```

System ID - TD-100

Run status - Run complete (this indicates a successful run)

Date of run

Time of run completion

Probe identification number

Operator identification number

TD-5 Lot #

Diagnostic status - System Test: Passed

Remarks:

In the remarks field you may add your own handwritten remarks for your records. For example, you may want to add a patient number or some other identification number.

## 4.2 Aborted/Invalid Cycle Printout

Should the system determine an error or incorrect procedure, such as an attempt to remove the probe during disinfection the system will abort the cycle. The system will print an aborted run report (as shown below) indicating a failed run.

```
TD-100
RUN ABORTED
*****
SYSTEM FAIL
WARNING
PROBE NOT DISINFECTED
*****

01/11/06
11:15:51A

PROBE #: 0
OPERATOR #: 0
TD-5 LOT #: 000

SYSTEM TEST: FAILED
REMARKS:

____
____
____
____
____
____
____
____

DIAG V1018:
01/11/06
11:15:51A
MODE: 021 000 004 000
A IN: 253 011 171 193 223
INIT: 253 011 000 171 194
TIME: 001 101 000 000

HIST: 000 045 041
M0: 001 020
M3: 004 000 000
M4: 006 004 008
M5: 003 000 000
M6: 024 000 000 002 006
M7: 000 000
```

The report at the left shows that an invalid run has occurred, “RUN ABORTED” indicates a failed run. The TEE probe is not disinfected.



**WARNING! THE TEE PROBE IS NOT DISINFECTED WHENEVER AN INVALID RUN OCCURS. THE PROBE IS CONSIDERED CONTAMINATED UNTIL A VALID RUN WITH A VALID RUN PRINTOUT OCCURS.**



## 4.3 Troubleshooting

For assistance with troubleshooting, please contact:

CS Medical  
3300 Lake Woodard Dr.  
Raleigh, NC 27604 USA

Phone: (877) 255-9472

Phone: (919) 255-9472

Fax: (919) 212-8050

Email: [service@csmedicalllc.com](mailto:service@csmedicalllc.com)

When you call for assistance, please have a copy of your invalid run printout.



**WARNING! REPAIRS BY UNAUTHORIZED PERSONEL SHOULD NOT BE ATTEMPTED AND WILL VOID THE MANUFACTURERS WARRANTY. UNAUTHORIZED REPAIRS MAY ALSO DAMAGE OR MAKE UNSAFE THE PRODUCT AND / OR ITS COMPONENTS.**

## 4.4 Water Supply Problem

If water pressure to TD-100™ drops below 10 PSI, the display will read LOW WATER ERR to indicate the water pressure is low. The unit will also audibly beep until the problem is corrected.

Check water supply

- Is the inlet hose connected?
- Is the supply line kinked?
- Is the pressure gauge at 20-25 PSI static?

## 4.5 Diagnostic Failure

During operation, the TD-100™ is constantly checking the valid function of the unit. If a fault is detected, the TD-100™ will shut down and an error report will print.

The diagnostic information at the bottom of the error report can be used to help diagnose the problem that occurred.

**NOTE: Should an error occur, please have a copy of the error printout available when you call customer service.**

In most cases, there are very few operator corrections that can be made when a diagnostic error has occurred. Any operator corrections that can be made are found in the following chart.

## 4.6 Troubleshooting

Problem/Message	Possible Cause	Corrective Action
Unit will not start / No display	Power switch off	Verify power switch is turned on
	Power cord loose or disconnected	Verify that power cord supplied with unit is firmly seated in unit and wall receptacle
	Tripped GFI outlet	Verify that power is available at outlet
	Blown fuse	Contact CS Medical
Display not legible	Display or controller problem	Contact CS Medical
Unusual or chattering noise	Controller problem	Contact CS Medical
Low water pressure message	Water supply line not connected	Verify that water line is connected per installation instructions
	Water supply line obstructed	Verify that water line is not pinched or kinked
	Water supply valve off	Verify that water supply valves are turned on
	Improperly adjusted regulator	Adjust to specifications
	Faulty pressure sensor	Contact CS Medical
Pump runs continuously when turned on	Drain not connected	Verify that drain line is connected per installation instructions
	Drain line obstructed	Verify that drain line is not pinched or kinked
	Dry unit	Pour one pint of water into reservoir. Water should pump immediately to drain and stop.
	Faulty valve or pressure sensor	Contact CS Medical
	Drain hose may have been shortened	Contact CS Medical
Liquid left in tube after run	Drain not connected	Verify that drain line is connected per installation instructions
	Drain line obstructed	Verify that drain line is not pinched or kinked
	Plugged strainer	Clean strainer per instructions.
	Faulty pump, valve or sensor	Contact CS Medical
Probe does not slide easily into tube	Incorrect probe	Verify that probe is approved for use in the TD-100™
	Obstruction in tube	Clean strainer per instructions.
Cycle will not progress beyond "close door" prompt	Door misaligned or obstructed	Visually inspect for obstructions and that door seats properly.
	Faulty switch	Contact CS Medical
Prompt to remove probe when no probe in place	Probe switch stuck	Press and release switch at rear of probe holder. Verify that switch moves freely.
	Faulty switch	Contact CS Medical
Prompt to remove TD-5™ when no TD-5™ in place	TD-5™ switch stuck	Press and release switch at rear of TD-5™ reservoir. Verify that switch moves freely.
	Faulty switch	Contact CS Medical

<b>Problem/Message</b>	<b>Possible Cause</b>	<b>Corrective Action</b>
No printout when run completes	No printer paper	Replace printer paper per instructions
	Printer paper jammed	Open printer door. Verify that paper runs through the slot in the printer door and is not jammed.
	Printer failure	Contact CS Medical
Water leakage observed around unit	Leaking water supply line or connection to unit	With water supply turned off, verify that quick disconnects are properly seated and hose clamps are tight.
	Disconnected or broken drain line	Verify that drain line is properly seated into correct port. Contact CS Medical for replacement if damaged.
	System overfilling during rinse	Contact CS Medical
	Loose or missing cleanout plug	Replace or tighten cleanout plug.
	Internal failure	Contact CS Medical
Disinfectant leakage	Leaking disinfectant bottle	Do not use leaking bottles. Clean up per instructions.
	Disconnected or broken drain line	Verify that drain line is properly seated into correct port. Contact CS Medical for replacement if damaged.
	Loose or missing cleanout plug	Replace or tighten cleanout plug.
	Internal failure	Contact CS Medical
Disinfectant overflow from reservoir	System not empty before starting cycle	See "Liquid left in tube after run"
Disinfectant odor	Air filter expired	Replace air filter per instructions
	Disinfectant leakage or spillage	Clean up per instructions.
	Doors not closed	Close doors
	Power interrupted during run	Restart unit.
	Fan failure	Contact CS Medical
Level error	Float switch stuck	Verify that float switches at rear of reservoir move freely up and down.
	Faulty pressure sensor	Contact CS Medical
Heater error	Low disinfectant	Verify that disinfectant is not leaking. If not then restart unit with full disinfectant bottle.
	Plugged strainer	Clean strainer per instructions.
	Faulty heater or thermal cutout	Contact CS Medical
	Faulty pump or valve	Contact CS Medical
Heater fail	Loss of heater control	Contact CS Medical

<b>Problem/Message</b>	<b>Possible Cause</b>	<b>Corrective Action</b>
Pump error	Plugged strainer	Clean strainer per instructions.
	Faulty pump, valve or sensor	Contact CS Medical
Pump fail	Loss of pump control	Contact CS Medical
Valve error	Drain not connected	Verify that drain line is connected per installation instructions
	Drain line obstructed	Verify that drain line is not pinched or kinked
	Faulty pump, valve or sensor	Contact CS Medical
Purge error	Plugged strainer	Clean strainer per instructions.
	Drain not connected	Verify that drain line is connected per installation instructions
	Drain line obstructed	Verify that drain line is not pinched or kinked
	Faulty pump, valve or sensor	Contact CS Medical
Fill error	Loose or missing cleanout plug	Replace or tighten cleanout plug.
	Started with disinfectant TD-5™ not full	Restart unit with full disinfectant bottle.
	Water supply line not connected	Verify that water line is connected per installation instructions
	Water supply line obstructed	Verify that water line is not pinched or kinked
	Faulty regulator, valve or sensor	Contact CS Medical
Overflow error	System not empty before beginning cycle	See "Liquid left in tube after run"
	Float switch stuck	Verify that float switches at rear of reservoir move freely up and down.
	Faulty valve or sensor	Contact CS Medical

## **SECTION 5: Maintenance**

### **5.1 Daily Checks and Daily Cleaning**

#### **Daily Checks**

Check all exterior plumbing connections on a daily basis and inspect for any sign of damage or leakage. Check the connections in the following order:

- Check the incoming water supply. Make sure the connection is not leaking and that the hose from the connection is not leaking or kinked.
- Check that the drain connection is secure and no leaks or damage are visible.
- Looking at the base or bottom of the unit observe there are no signs of leakage or “spots” which might indicate prior leakage.

Should any damage or leaks be found, refer to the troubleshooting section of this manual to resolve the problem.

#### **Daily Cleaning**

Clean all external surfaces with a soft cloth dampened with 70% Isopropyl Alcohol. Open the probe door and the TD-5™ door to clean internal surfaces with a soft cloth dampened with 70% Isopropyl Alcohol.

## 5.2 Changing the Printer Paper

To change the printer paper, open the cover by pulling up on the colored tab, see figure 14. To avoid damage do not use excessive force.



Figure 14, Opening Printer Paper Cover.




This will expose the chamber that holds the thermal paper roll. Remove and discard the roll. Place a new roll into the chamber paying close attention to the direction of installation. The roll should be placed in so the paper is feed from the bottom of the roll.

Pull out an inch or two of paper and close the cover by applying equal amounts of pressure on each side until the cover locks in place, see figure 15. Tear away the excess paper and the printer is ready to use.



Figure 15, Closing Printer Paper Cover.

### 5.3 Changing the Air Filter (Part # CS-TDFGLT)

-  **CAUTION! THE AIR FILTER MUST BE CHANGED ACCORDING TO YOUR SERVICE AGREEMENT.**
-  **CAUTION! TURN OFF THE POWER BEFORE REPLACING THE AIR FILTER.**
-  **CAUTION! FILTER IS HEAVY.**



To change the filter:

1. Remove filter door.
2. Remove the filter.
3. Replace with new filter taking care not to damage the gasket on the filter as you insert it.
4. Ensure arrows are pointing upwards.
5. Replace cover. Unit is now ready to use.

Fig. 16



## 5.4 Setting or Changing the Date and Time.

When you first receive your TD-100™, you may need to set the date and time. Follow the steps below to enter the correct date and time. Do not press “ENTER” except at the end of the sequence below.

- 1) Turn on the main power to the TD-100™ by activating the main power switch.
- 2) Unit will briefly enter the purge mode and then display “PRESS START” on the top line of the display and “SETUP” on the bottom line.
- 3) Now press the “A” key to enter the setup mode. The TD-100™ will display a random time and date.
- 4) The date is displayed on the top line in MM/DD/YY format. The month should be blinking.
- 5) Enter a two digit number corresponding to the current month, e.g. 07 = July.
- 6) The TD-100™ will now automatically step to the day field. Enter a two-digit number corresponding to the current day.
- 7) The TD-100™ will step automatically to the year field. Enter the current year, e.g. 03 = 2003
- 8) The TD-100™ now steps to the time field. Time is displayed in HH/MM/SS AM or PM.
- 9) Enter a two-digit number corresponding to the current hour. The TD-100™ now automatically steps to minutes.
- 10) Enter a two-digit number corresponding to the current minute, e.g. 32.
- 11) The TD-100™ now steps to seconds. Enter a two-digit number corresponding to current second.
- 12) The TD-100™ now automatically assumes it is AM and displays “A” To change to PM press “1” on the keypad.
- 13) Now press, “ENTER” to accept your changes. TD-100™ displays “PRESS START” and “SETUP.”

NOTE: If you press “ENTER” by mistake before you are finished, select “SETUP” and go back to (1) above.

## SECTION 6: Warning and Caution Statement Summary

### 6.1 Warning Statements

Throughout this manual WARNING and CAUTION statements have been used to alert the operator to important information regarding operator safety or the correct use of this device to achieve high-level disinfection. Please be familiar with these statements before attempting to use the TD-100™ TEE probe disinfectant.

A summary of these statements is as follows:



#### Warning Statements



**These instructions should be placed in a visible location.**



**Do not use the TD-100™ until it has been properly installed.**



**Make sure the TD-100™ is connected to a hospital grade GFI electrical outlet only.**



**The TD-100™ is not designed to be used in explosion endangered areas.**



**The TD-100™ is not protected against water spillage or spray. Take precautions in the use of liquid during cleaning or insertion of devices.**



**Do not attempt to manually open a TD-5™ disinfectant container**



**Avoid contact with skin and avoid breathing vapors of the TD-5™ disinfectant.**



**The TD-100™ may only be used with Transesophageal Ultrasound probes. Do not use any other device in this system.**



**Do not use any other disinfectant other than the TD-5™ disinfectant.**



**Do not attempt to refill the TD-5™ with any type of disinfectant.**



**Follow maintenance procedures in the Operators Manual. Failure to follow maintenance procedures may compromise the use of this device.**



**Always verify that the disinfectant container is empty at the completion of a cycle.**



**Repairs by unauthorized individuals should not be attempted and may result in damage or malfunction. Non-operator serviceable components should be serviced by an authorized representative of CS Medical.**

## SECTION 7: Spare Parts

### 7.1 Spare Parts List

Part #	Description
200293	Power Cord*
200496	Water Supply Hose*
200497	Drain Hose*
200499	Drain Filter*
CS-TDFGLT	Air Filter
200295	Printer Paper
200258	Fuse, 10A

\* See Installation Instructions for Replacement

## **SECTION 8: Contact Information**

### **8.1 Contact Information**

Manufactured by:

CS Medical LLC  
3300 Lake Woodard Drive  
Raleigh, NC 27604 USA

[www.csmedicalllc.com](http://www.csmedicalllc.com)

Phone: (877) 255-9472

Phone: (919) 255-9472

Fax: (919) 212-8050